



Lightning-Link™ Biotin Conjugation Kit (Type B*)

***Optimized for assays in which conjugate is captured by a streptavidin plate**

715-0010 3 x Ab labelings (each 100-200mg scale)
715-0015 1 x Ab labeling (1-2 mg scale)

Release 001; revised June 2007


Technical bulletin 1097

1. INTRODUCTION

The Lightning-Link™ conjugation kit allows biotinylation to be set up in *seconds*, simply by adding a solution of the protein to be labeled to a lyophilised mixture containing a proprietary activated biotin ligand (Figure 1).

By circumventing the desalting or dialysis steps that commonly interrupt traditional protein conjugation procedures, Lightning-Link™ technology can be used to label small quantities of protein with 100% recovery and no excessive dilution of the conjugate.

Fig. 1 Lightning-Link™ antibody conjugation

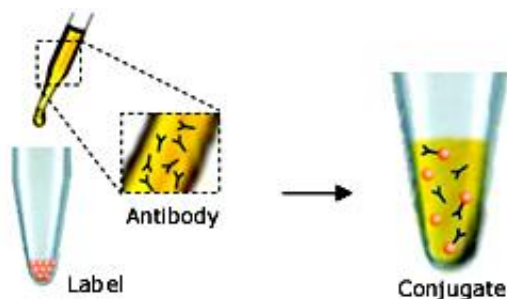

Upon dissolution of Lightning-Link™ mixture with a solution of the antibody (or other biomolecule to be labeled) proprietary chemicals in the mixture become activated. This results in coupling of the antibody to the biotin, which has an extended spacer, in a gentle and controlled process at near-neutral pH. Lightning-Link™ makes it possible to biotinylate primary antibodies and other proteins with ease, using a simple, highly efficient process that has safeguards to prevent over-labeling of the biomolecule and that ensures 100% recovery even at small scale.

2. INSTRUCTIONS

2.1 Storage and components

The kit is shipped at ambient temperature in a tamper-evident polypropylene container. Store at -20°C upon receipt.

Kit contents:

1 or 3 glass vial(s) of Lightning-Link™ mix

1 vial of LL-Modifier reagent

1 vial of LL-Quencher reagent

2.2 Considerations before use

2.2.1 Sample buffer

Ideally, the antibody to be labeled should be in 10-50mM amine-free buffer pH range 6.5 to 8.5. However, many buffers outside these limits of concentration and pH can be accommodated. Modest concentrations of Tris buffer are also tolerated. Appendix 1 gives further guidance on buffers and compatible additives.

2.2.2 Amount and volume of antibody

The recommended amount of antibody to be used for biotinylation reactions is up to 200µg for 715-0010 and up to 2mg for 715-0015. The volume of the antibody sample, ideally, should be in the range 40-100µl (715-0010), and 400-1000ul (715-0015). Antibody concentrations 1mg/ml or greater generally give the best results. The LL-biotin kit contains a stable form of biotin and is much more efficient than standard biotin kits, thus no purification of the conjugate is required (see also FAQ 2).

2.3 Setting up conjugation reactions

2.3.1. Before you add antibody to the Lightning-Link™ mix, add 1µl of LL-Modifier reagent for each 10µl of antibody to be labeled. Mix gently.

2.3.2. Remove the screw cap from the vial of Lightning-Link™ mix and pipette the antibody sample (with added LL-modifier) directly onto the lyophilised material. Resuspend *gently* by withdrawing and re-dispensing the liquid once or twice using a pipette.

2.3.3. Place the cap back on the vial and leave the vial standing for 3 hours at room temperature (20-25°C). Alternatively, and sometimes more conveniently, conjugations can be set up and left overnight, as the longer incubation time has no negative effect on the conjugate.

2.3.4. After incubating for 3 hours (or more), add 1µl of LL-quencher reagent for every 10µl of antibody used. The conjugate can be used after 30 minutes.

2.4 Storage of conjugates

For any new conjugate, initial storage at 4°C is recommended. A preservative may be desirable for long-term storage. Other storage conditions (e.g. frozen at -70°C or stored at -20°C with 50% glycerol) may also be satisfactory. The best conditions for any particular conjugate must be determined by experimentation.

Appendix 1. Compatibility of buffers and buffer additives.

Amine-free buffers, including MES, MOPS, HEPES and phosphate are compatible if they are in the concentration range 10-50mM and have pH values in the range 6.5-8.5, as the addition of LL-Modifier provides the conditions necessary for efficient conjugation. Common non-buffering salts (e.g. sodium chloride), chelating agents (e.g. EDTA), and sugars may be present, as they have no effect on conjugation efficiency. Azide (0.02-0.1%) has little or no effect.

If the amine-free buffer is relatively concentrated and outside the pH range 6.6-8.5 you may need to add more LL-modifier for each 10µl of antibody. Excess LL-Modifier is provided so that you can check the pH of the buffer after the addition of the modifier. Ideally the pH should be around 7.3-7.6, though efficient conjugation occurs anywhere between pH 6.8 and 7.8. Avoid

buffer components that are nucleophilic, as these may react with Lightning-Link™ chemicals. Primary amines (e.g. amino acids, ethanolamine) and thiols (e.g. mercaptoethanol, DTT) fall within this class. (Note: Tris has little effect on conjugation efficiency as long as the concentration is 20mM or less).

Q1. What functional groups do I need on my protein?

Lightning-Link™ requires amine groups on the molecule to be labeled. Most proteins have lysine and/or alpha-amino groups. All antibodies will have multiple amine functions.

Q2. Do I need to purify the conjugate?

No. The chemicals used in Lightning-Link™ are deactivated by the quencher, and the by-products are benign. Moreover, the conjugation efficiency is very high and, unlike other biotin reagents, LL-biotin does not hydrolyze in solution or during storage, thus it is not necessary to use a large molar excess or to purify the final conjugate.

Q3. Can non-antibody molecules be labeled?

Yes. While labeling of antibodies is a major application, the only requirement is that the protein to be labeled has amine functionality.

Limited use license

Innova Biosciences' Lightning-Link™ conjugation kits are offered for research purposes alone, and are not intended for human, therapeutic or diagnostic use. The purchase of this conjugation kit conveys to the buyer (whether the buyer is a not-for-profit, academic or for-profit entity) the non-transferable right to use the amount of product purchased and the components of the product for in-house research. The buyer shall not sell or otherwise transfer this product, its components, or materials prepared therefrom to any third party. The buyer shall not use this product or its components for commercial purposes. For the avoidance of doubt, 'commercial purposes' means any activity by a party for consideration and includes, without limitation, use of the product or its components (i) in the manufacturing of conjugated materials (e.g. labeled antibodies), (ii) to provide a service, information or data, (iii) for therapeutic, diagnostic or prophylactic purposes, or (iv) for repackaging/resale, whether or not such product or its components are resold for use in research. The use of this product by the buyer constitutes agreement with the terms of this limited use label license for Lightning-Link™ products.

For information on purchasing a license for commercial applications contact Innova Biosciences Ltd, Business Development Office, Babraham Hall, Babraham, Cambridge, UK, CB22 3AT. Tel +44(0)1223 496170; Fax +44(0)1223 496172.